


1

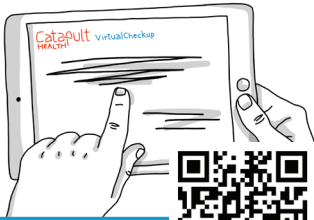

WATCH THE VIDEO

- Before you begin, please watch our 4-minute “how to” video at VirtualCheckup.com/video
- You may also scan this QR Code with your phone’s camera to access the video.
- Drink plenty of water before you begin.


2



BEGIN CONSENT FORM

- The information you include on this form is confidential.
- All fields are required.
- The more accurate and honest you are, the more we can help you!
- Weigh yourself, record your height and weight.
- Answer the question about fasting.


3

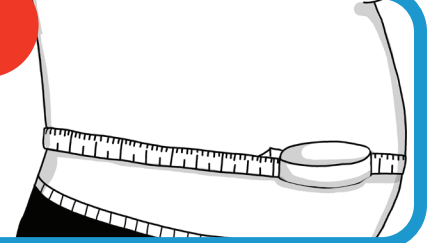

MEASURE YOUR ABDOMEN

- Locate measuring tape in kit.
- Measure the distance around your abdomen, right above your hip bone and around your belly.
- Your abdominal measurement is not the same as your pant size.
- Abdominal circumference is often greater than pant size.
- Record on consent form.


4

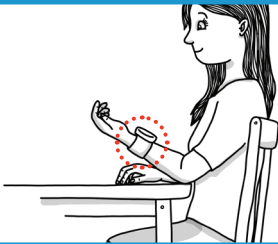

CHECK YOUR BLOOD PRESSURE

- Sit in a chair pulled up to a table.
- Place your feet flat on the floor.
- Back straight, elbow on the table.
- Install batteries in blood pressure monitor, put on your left wrist.
- Screen should be visible with your palm up (opposite to wearing a watch).
- Relax wrist and hand, raise wrist to heart level.
- Take blood pressure reading, record on consent form, then wait 3 minutes.
- Take second reading and record on consent form. Also record heart (pulse) rate from monitor.


5

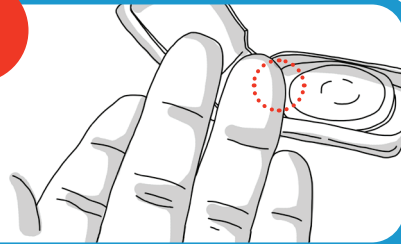

GET EVERYTHING READY

- Open the foil pouch and remove the green plastic collection device.
- Place collection device (with collection side up) on a flat and level surface.
- Do NOT remove the clear plastic that covers the white collection site.
- Rinse hands in warm water.
- With palm facing up, select ring finger for your collection. Puncture site should be on the **side of the fingertip** for ease of collection. Do NOT use the center of the fingertip.


6

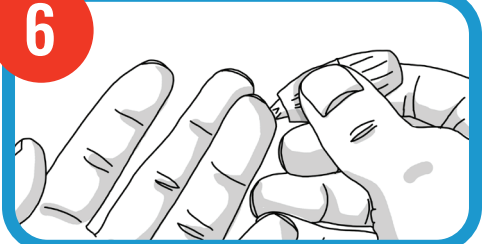

THE FINGER STICK MADE EASY

- Remove ring(s) from ring finger.
- Clean collection site with alcohol pad.
- Stimulate blood flow by letting your hand hang at your side and shaking back and forth for 15-20 seconds.
- Twist off cap of lancet provided. Place hand, palm side up, on a hard surface.
- Position lancet on right side of finger as shown. Press lancet firmly into side of fingertip until it “clicks.”
- Use gauze to wipe off and discard first drop of blood.


7

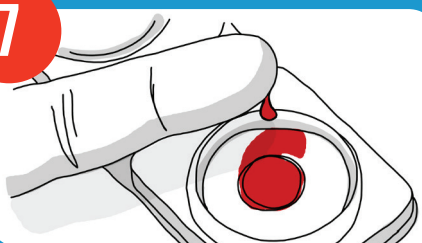

4 DROPS IS ALL IT TAKES

- Squeeze above the knuckle of the finger to get the collection started.
- Blood drops should flow freely.
- If blood drop does not form, wipe site with gauze pad and gently massage finger from the palm to fingertip.
- **Four (4) large drops are required** to fall through the center hole of the clear cover onto the application surface.
- Apply drops continuously. Do NOT STOP or wait between drops.
- When finished, hold gauze on finger 1 minute, remove and apply bandage.

8

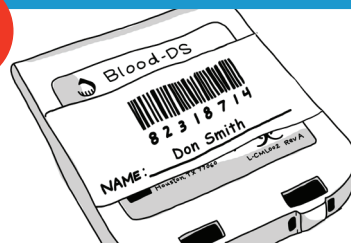

CLOSE AND LABEL THE SAMPLE

- After 4 blood drops are absorbed, wait 3 minutes then close the device firmly by folding over the top.
- Make sure the latches are secure.
- Remove one barcode label from the consent form, print your name on it, affix to the outside of the collection device.
- The collection device containing your blood sample can be stored briefly at room temperature (65°-72°) before mailing to the laboratory.

9



MAIL YOUR SAMPLE & FORM

- Verify all information is correct and complete on the consent form. Sign and date.
- Remove a second barcode label and affix to the Participant Card located in the box tray. Retain this card.
- **Place your blood collection device and your signed consent form in the prepaid return envelope provided.**
- Your blood sample and consent form must be mailed (US MAIL) within 24 hours of collection.

What's Included in your Home Kit?

All fields on your Test Requisition and Consent Form **must be completed**. Form must be returned with your blood sample.



DO NOT OPEN foil pouch until you're ready to begin collecting your blood sample. The collection device inside the pouch **must be used within 30 minutes** after opening foil pouch.

WARNING: Important Collection Information

Improper collection of blood specimens may affect laboratory results. Carefully read and follow the collection instructions. Persons with hemophilia or anti-coagulation therapy should consult a licensed physician or healthcare professional before using this kit.


Before you begin, please watch our video at VirtualCheckup.com/video

You may also scan this QR Code with your phone's camera to access the video.



Catapult **VirtualCheckup**[™]
HEALTH

200 NE Missouri Road, Suite 304 Lee's Summit, MO 64086 | coremedicalabs.com

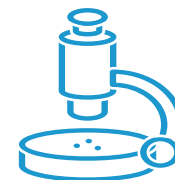

Welcome to
Catapult **VirtualCheckup**[™]
HEALTH

Convenient Self-Collection Kit